

If there is a fault with the appliance, inform: _____

Model ID / _____

Inventory number: _____

Quick Start Guide for everyday use

Information about using this Quick Start Guide:

- This Quick Start Guide is not a substitute for the full operating instructions nor is it a substitute for appliance training.
- This Quick Start Guide provides support for everyday use and handling the appliance.
- Read the full operating instructions and the safety notes contained therein.
- This Quick Start Guide is for: SRFvh 40../55..
- Scanning the QR code will open the full operating instructions:



Intended use

This laboratory refrigerator is suitable for professional storage of products at temperatures between 3 °C and 16 °C.

Foreseeable incorrect use

Do not use the appliance for the following applications:

- Storage and refrigeration of:
 - Chemically unstable, flammable or corrosive substances
 - Blood, plasma or other body fluids for the purpose of infusion, application or introduction in the human body
- Use in potentially explosive atmospheres
- Use outdoors or in areas exposed to dampness and splashing water

When stocking, observe the following points:

- ☐ Position grid shelves according to the height required.
- ☐ Observe maximum load.
- ☐ When the storage temperature has been reached, load the appliance (compliance with cold chain).
- ☐ Keep the ventilation slots of the recirculation fan free in the interior.
- ☐ Produce for cooling must not come into contact with the evaporator on the rear wall.
- ☐ Keep liquids in closed containers.
- ☐ Leave space when storing refrigerated goods to ensure good air circulation.

Status display

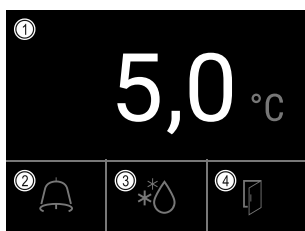


Fig. 155

(1) Actual temperature
(2) Alarm status

(3) Defrosting status
(4) Door opening status

Navigation with the touch display

Navigation	Description
	Swipe right or left Navigates within the menu.
	Press briefly Activates/deactivates function. Confirms selection. Opens submenu.
	Long press (3 seconds) Activates/deactivates function.
	Brief press of Back symbol Jumps back one level.

Displaying temperature recording

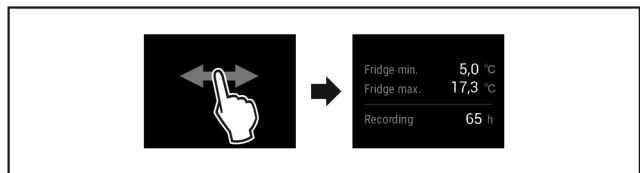


Fig. 156

► Complete steps (see Fig. 5562) .

Resetting temperature recording

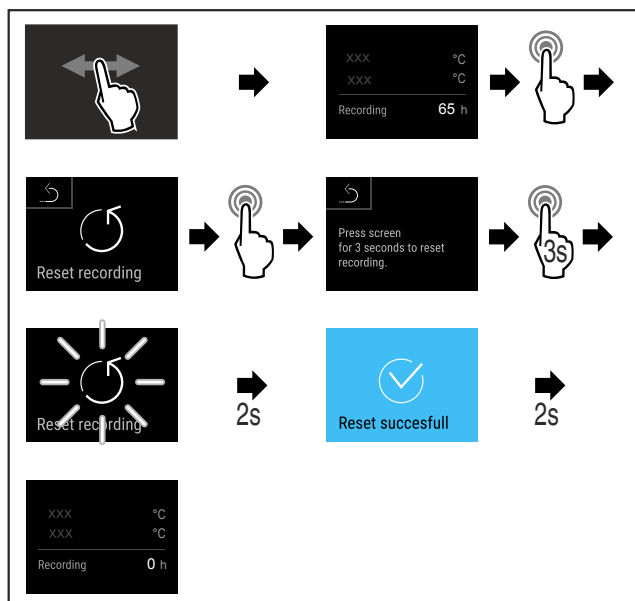


Fig. 157

► Complete steps (see Fig. 5563) .

Quick Start Guide for everyday use

Unlocking appliance

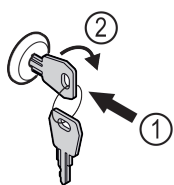


Fig. 158

- Insert the key. Fig. 158 (1)
- Turn the key clockwise by 90°. Fig. 158 (2)
- ▷ The appliance is unlocked.

Displaying alarm log

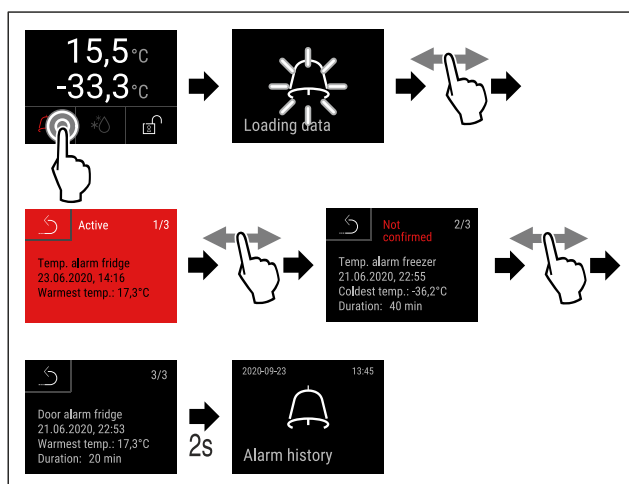


Fig. 159 Example shows an active alarm, an unconfirmed alarm, and a confirmed alarm.

- Complete steps (see Fig. 5565) .
- ▷ Alarm log shows active, confirmed or unconfirmed alarms (maximum 10 alarms).
- ▷ To delete the alarm log: See operating instructions.

Cleaning appliance

- Empty the appliance.
- Pull out the mains plug.
- Observe the information in the operating instructions.
- Dry the appliance and equipment thoroughly after cleaning.
- Connect the appliance.
- Switch on the appliance.
- Load the appliance once the storage temperature is reached.



What should I do if there is a fault?

A message appears in the display if there is a fault. Messages have different meanings:

Symbol colour in the display	Meaning	Remedy
YELLOW	Reminder Reminder of general procedures. You can carry out these procedures to eliminate the message.	Consult the operating instructions and perform the troubleshooting steps.

Symbol colour in the display	Meaning	Remedy
RED	Warning Appears when there are malfunctions. You can eliminate simple malfunctions yourself.	Consult the operating instructions and perform the troubleshooting steps.
If you cannot remedy the fault yourself, contact a competent member of staff or your nearest customer service outlet.		











Cut out the page



Quick Start Guide for everyday use

Functions

Overview of functions

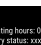	Switch appliance on and off	 Door alarm	Door alarm
	Reset recording	 Light alarm	Light alarm
	Manually starts automatic defrosting	 Alarm simulation	Alarm simulation
	Settings Access to this function can be blocked using a PIN code. (see 5.3.6 Settings menu access protection)	 Alarm relay	Alarm relay
	Info	 Alarm history	Alarm log
	Temperature setting	 Calibration	Sensor calibration
	Temperature alarm	 Maintenance timer	Maintenance interval
	Activate/deactivate settings lock	 Light	Lighting
	Alarms	 Language	Language
	Button sounds	 Date / time	Date/time
	Door light	 Temperature unit	Temperature unit
	Defrost information	 Display brightness	Screen brightness
	Sensor Info	 Download data	Data download
	DemoMode	 LAN	LAN *
	Factory reset	 WiFi	WiFi *
		 System status	System status
		 Software	Software
		 Settings code	Settings lock PIN code change