

Electronic Stethoscope

Powerful Capability. Incredible Simplicity.
Redefining what a stethoscope can do for you.

Model 3100

- State-of-the-art sound sensor (offers the sound of a high-end cardiology stethoscope)
- Ambient noise reduction reduces approx 85% of distracting environmental noise
- Up to 24 X sound amplification
- LCD interface to keep you informed about current settings and battery life
- Worry-free power system that includes multiple reminders
- Versatile chestpiece can accommodate infant, paediatric or adult patients
- Two year warranty



Model 3100

- Black 3M3100BK
- Navy Blue 3M3100NB
- Burgundy 3M3100BU
- Plum 3M3100PL

\$495



Model 3200

Using wireless bluetooth technology, the 3200 model is a highly capable yet user friendly stethoscope. It includes all the features of the Model 3100 plus more:

- Bluetooth® technology wirelessly transmits heart and lung sounds to your PC for further analysis or transfer to patient records
- Built in recording feature records up to 12 sounds for playback
- Provides three filter modes for optimal heart and lung auscultation - Bell 20–200 Hz, Diaphragm 100–500 Hz, Extended Range 50–500 Hz
- Includes free heart and lung sound visualisation software useful for detecting heart murmurs
- Two year warranty

Model 3200

- Black 3M3200BK
- Navy Blue 3M3200NB
- Burgundy 3M3200BU

\$619

Littmann Stethoscope spare parts

Diaphragm and Rims			
a	Cardiology III and Classic II	Black	3M36556
	Cardiology III paediatric side	Black	3M36572
Ear Tips			
b	Large, soft, snap on	Black	3M37811
	Small, soft, snap on	Black	3M37809
Rims			
c	Cardiology III large side and Classic II	Black	3M36548
	Classic II Paediatric and Cardiology III small side	Black	3M36550
Non-Chill Bell Sleeves			
d	Cardiology II and Cardiology III	Black	3M36570
	Classic II	Black	3M36544



The full range of Littmann spares is available

Littmann Classic II Paediatric

Designed and sized to provide Littmann quality acoustics for children

- Double-sided stainless steel chestpiece
- 2.5cm traditional paediatric bell
- Single lumen binaural
- Three year warranty

- Black 3M2113
- Red 3M2113R
- Carribean Blue 3M2119
- Raspberry 3M2122

\$119



Littmann Classic II S.E.

The industry standard

- High acoustic sensitivity
- Traditional bell combined with a tunable diaphragm allowing the user to alternate between low- and high-frequency sounds
- Ideal for daily use by General Practitioners, Nurses
- Three year warranty

- Black 3M2201
- Burgundy 3M2211
- Navy Blue 3M2205
- Purple 3M2209
- Orange 3M2812
- Carribean Blue 3M2206
- Raspberry 3M2210
- Ceil Blue 3M2813
- Bubblegum Pink 3M2817
- Hunter Green 3M2208
- Seafoam Green 3M2814
- Pine Green 3M2818
- Ocean Blue 3M2819
- Grey 3M2203

\$99



Littmann Cardiology III

Outstanding acoustic performance and exceptional versatility

The Littmann Cardiology III stethoscope offers two tunable diaphragms which allow the user to conveniently alternate between low- and high-frequency sounds without turning over the chestpiece. The large side can be used for adult patients, while the small side is especially useful for paediatric or thin patients, around bandages and for carotid assessment.

The paediatric side of the chestpiece easily converts to a traditional bell by simply replacing the diaphragm with the non-chill bell sleeve included with each stethoscope.

- Five year warranty

- Black 3M3128
- Navy Blue 3M3130
- Burgundy 3M3129
- Plum 3M3135

\$239



Littmann Master Cardiology

The Master Cardiology meets the standards of healthcare professionals who demand very high acoustic response. It offers a handcrafted stainless steel chestpiece with the portable convenience of a single-sided stethoscope. The tunable diaphragm conveniently alternates between low- and high-frequency sounds with a simple pressure change on the chestpiece.

- Superior performance
- Tunable diaphragm alternates between low- and high-frequency sounds
- Seven year warranty

- Black 3M2160
- Navy Blue 3M2164
- Burgundy 3M2163
- Plum 3M2167

\$269





Classic stethoscope

Chestpiece and binaural made from durable surgical grade stainless steel



- Black [BYDI417BK](#)
- Navy Blue [BYDI417NABL](#)
- Hunter Green [BYDI417HUNGN](#)
- Grey [BYDI417GY](#)
- Light Blue [BYDI417LIBL](#)
- Burgundy [BYDI417BUR](#)

\$64.95



Cardiology stethoscope

- Extra deep bell complete with non-chill ring for patient comfort
- Super sensitive diaphragm for low frequency sounds.



- Black [BYDI421BK](#)
- Green [BYDI421GN](#)
- Grey [BYDI421GY](#)
- Burgundy [BYDI421BUR](#)

\$89



Sprague rappaport

- Unlimited flexibility and superior acoustics
- Five interchangeable chestpiece fittings
- 22" heavy walled double tube configuration for superior acoustic separation

- Black [BYDI418BK](#)
- Grey [BYDI418GY](#)
- Burgundy [BYDI418BUR](#)
- Blue [BYDI418BL](#)


\$34.95



- Black [BYDI415BK](#)
- Blue [BYDI415BL](#)
- Green [BYDI415GN](#)
- Grey [BYDI415GY](#)
- Pink [BYDI415PK](#)
- Purple [BYDI415PUR](#)
- Red [BYDI415RD](#)
- Yellow [BYDI415YW](#)

Single Head Nurse Stethoscope

\$12.95

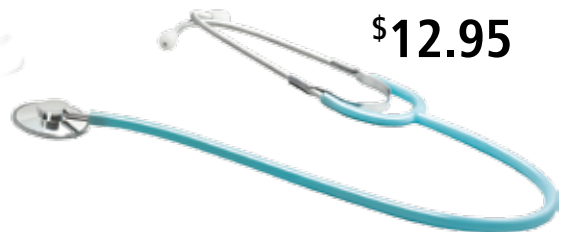

Dual Head Stethoscope



\$22.95

- Black [BYDI416BK](#)
- Blue [BYDI416BL](#)
- Green [BYDI416GN](#)
- Grey [BYDI416GY](#)
- Pink [BYDI416PK](#)
- Purple [BYDI416PUR](#)
- Red [BYDI416RD](#)
- Yellow [BYDI416YW](#)



Dual Head Stethoscope

\$12.95

[AXP10111511](#)



Toonscope

- Single head stethoscope with distinctive molded resin animal face chestpiece
- 6 interchangeable animal face snap-ons included
- Lightweight brushed aluminium binaurals
- Includes spare ear tips
- Comes in 4 colours

- Light Blue **BYDI425LB**
- Yellow **BYDI425YW**
- Lime **BYDI425LM**
- Pink **BYDI425PK**

\$49



Pinnard Foetal Stethoscope



- Aluminium **BYDS426** **\$7.95**
- Plastic **BYDS420** **\$6.95**

Training Stethoscope



BYDI565 **\$32.95**

HEINE Gamma Stethoscopes

HEINE GAMMA 3.1® Pulse Stethoscope

For fine acoustic performance

- Stethoscope with solid flat chest piece. The chest piece can be easily inserted under the sphygmomanometer cuff
- High-quality membrane, 44mm
- Single-tube construction with excellent sound transmission

HEM00009941

\$78.95



HEINE GAMMA 3.2® Acoustic Stethoscope

For excellent acoustic performance

- Similar features as GAMMA 3.1, but with a double chest piece
- Combined double chest piece with two different frequency ranges to choose from
- High-quality membrane for excellent acoustic performance and sound transmission
- Cone diameter 30mm

HEM00009942

\$103.95



HEINE GAMMA 3.3® Acoustic Stethoscope

For the examination of children

- Similar features as GAMMA 3.2. but with a reduced contact area of 31mm
- Ideal for examining children

HEM00009943

\$103.95



HEINE GAMMA C3® Cardio Stethoscope

Specially-developed for cardiology

- Exceptional acoustic performance, particularly in the cardiac frequency range
- Combined double chest-piece with two different frequency ranges
- High-performance double membrane, 45mm diameter gives ideal sound transmission and perfect skin contact for exceptional sound reproduction
- Two-in-one tube system for enhanced performance

HEM00009944

\$289.95

