


header

Patient ID	1234567	Weight	80 kg	Height	175 cm
MRN	6789	Address	Second street, 200, Big City		
Last Name	Doe	Order	555667		
First Name	John	Print Date	08-Mar-23		
Birth Date	01-Jan-70	Age	53	Test Date	05-Aug-19
Gender	Male	Recorder model	ABPM Device	S/N	ENG356

Physician Comments

Physician can write here some comments

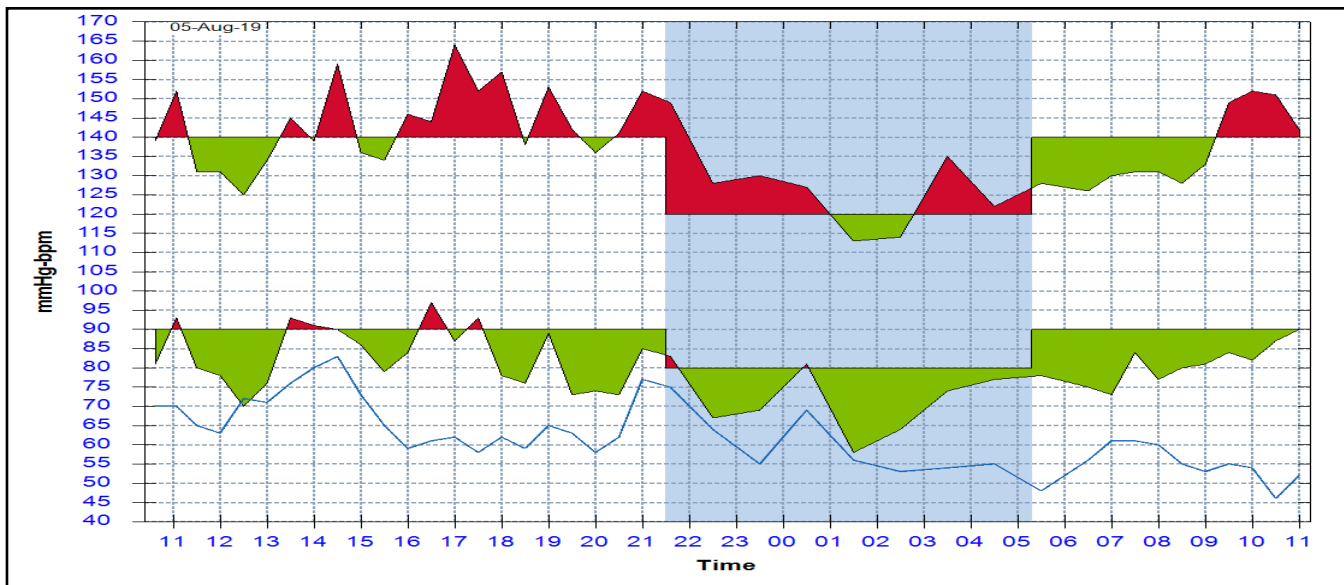

Brachial BP Results

Period	Time	Samples	Average SYS mmHg (+/-Std.Dev)	Average DIA mmHg (+/-Std.Dev)	Average HR BPM (+/-Std.Dev)	BP Load Sys (%)	BP Load Dia (%)
Overall	10:37-12:00 (25:23)	48	138 (+/-11.9)	80 (+/-8.4)	62 (+/-8.7)	54	22
Awake Period	05:18-21:30 (16:12)	40	141 (+/-10.5)	82 (+/-7.1)	63 (+/-8.9)	48	21
Asleep Period	21:30-05:18 (07:48)	9	127 (+/-11.6)	72 (+/-8.7)	60 (+/-8.2)	75	25
Dipping	SYS = 9.71% DIA = 13.01%						

Referring MD

Reporting MD

Signature

Date

Signature

Date

Patient ID	1234567	Weight	80	kg	Height	175	cm
MRN	6789	Address	Second street, 200, Big City				
Last Name	Doe	Order	555667				
First Name	John	Print Date:	08-Mar-23				
Birth Date	01-Jan-70	Age	53	Test Date:	05-Aug-19		
Gender	Male	Recorder model	ABPM Device	S/N	ENG356		

Medical History

From	To	Disease	Notes
07-Feb-23	08-Mar-23	Disease First	Some Notes 11111111

Medications


From	To	Trade Name	Active Agent	Dosage	Notes
07-Feb-23	08-Mar-23	Medication A	aaa	100	Some notes for the medication

header

Patient ID	1234567	Weight	80 kg	Height	175 cm
MRN	6789	Address	Second street, 200, Big City		
Last Name	Doe	Order	555667		
First Name	John	Print Date	08-Mar-23		
Birth Date	01-Jan-70	Age	53	Test Date	05-Aug-19
Gender	Male	Recorder model	ABPM Device	S/N	ENG356

Physician Comments

Physician can write here some comments

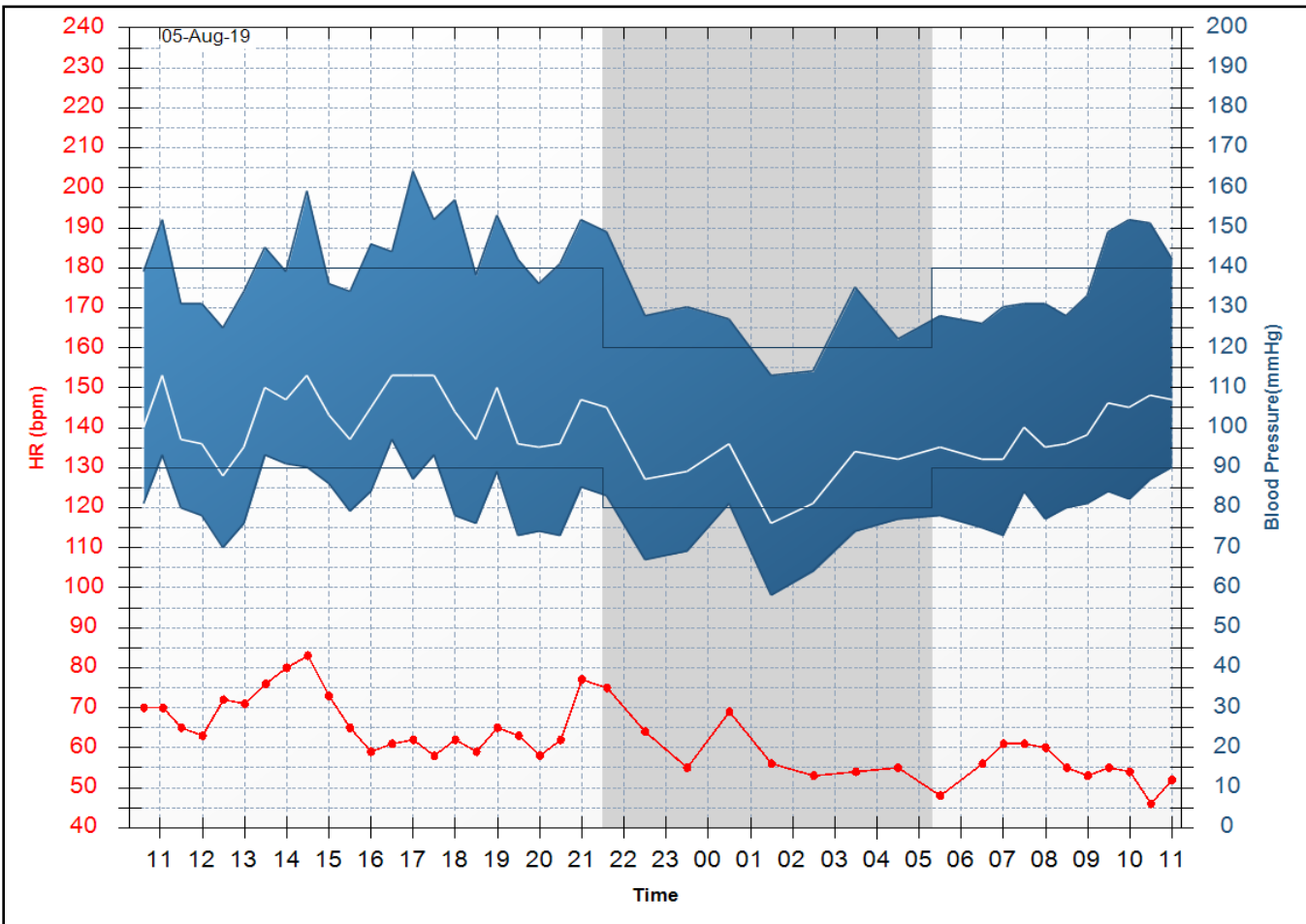


header

Patient ID	1234567	Weight	80 kg	Height	175 cm
MRN	6789	Address	Second street, 200, Big City		
Last Name	Doe	Order	555667		
First Name	John	Print Date	08-Mar-23		
Birth Date	01-Jan-70	Age	53	Test Date	05-Aug-19
Gender	Male	Recorder model	ABPM Device	S/N	ENG356

Physician Comments

Physician can write here some comments

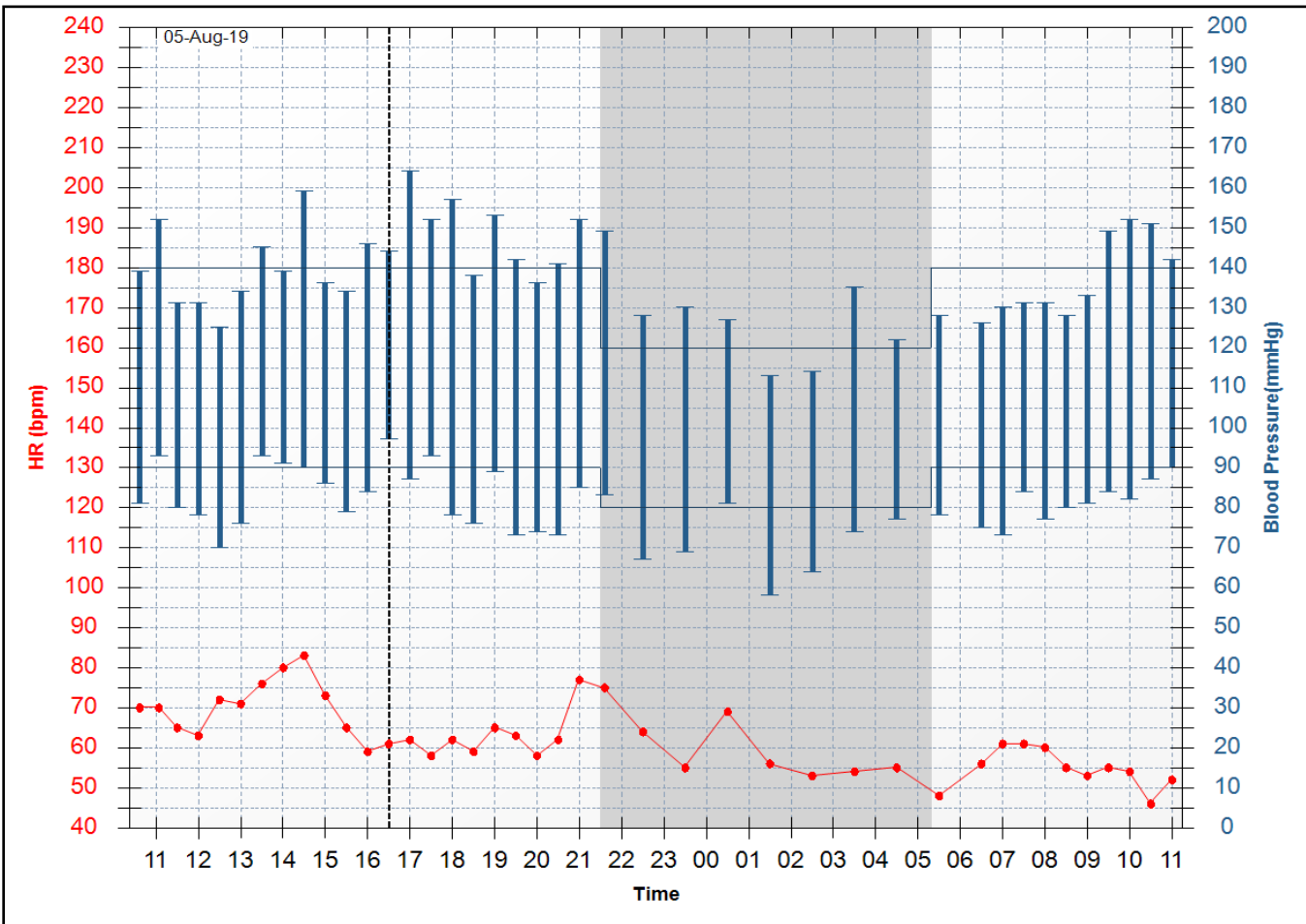

header

Patient ID	1234567	Weight	80	kg	Height	175	cm
MRN	6789	Address	Second street, 200, Big City				
Last Name	Doe	Order	555667				
First Name	John	Print Date	08-Mar-23				
Birth Date	01-Jan-70	Age	53	Test Date	05-Aug-19		
Gender	Male	Recorder model	ABPM Device	S/N	ENG356		

Physician Comments

Physician can write here some comments

Measurements

#	Test Date	SYS	MAP	DIA	HR	Code	Notes
1	05-Aug-19 10:37:00 AM	139	100	81	70	M	
2	05-Aug-19 11:00:03 AM	-	-	-	-	E 91	Service required (safety override is fitted or autozero is out of range).
3	05-Aug-19 11:00:05 AM	-	-	-	-	E 6	Release interval violation
4	05-Aug-19 11:04:07 AM	152	113	93	70	A	
5	05-Aug-19 11:30:03 AM	131	97	80	65	A	
6	05-Aug-19 12:00:03 PM	131	96	78	63	A	
7	05-Aug-19 12:30:03 PM	125	88	70	72	A	
8	05-Aug-19 1:00:03 PM	134	95	76	71	A	
9	05-Aug-19 1:30:03 PM	145	110	93	76	A	
10	05-Aug-19 2:00:03 PM	139	107	91	80	A	
11	05-Aug-19 2:30:03 PM	159	113	90	83	A	
12	05-Aug-19 3:00:03 PM	136	103	86	73	A	
13	05-Aug-19 3:30:03 PM	134	97	79	65	A	
14	05-Aug-19 4:00:03 PM	146	105	84	59	A	
15	05-Aug-19 4:30:03 PM	144	113	97	61	A	
16	05-Aug-19 5:00:03 PM	164	113	87	62	A	
17	05-Aug-19 5:30:03 PM	152	113	93	58	A	
18	05-Aug-19 6:00:03 PM	157	104	78	62	A	
19	05-Aug-19 6:30:03 PM	138	97	76	59	A	
20	05-Aug-19 7:00:03 PM	153	110	89	65	A	
21	05-Aug-19 7:30:03 PM	142	96	73	63	A	
22	05-Aug-19 8:00:03 PM	136	95	74	58	A	

header


23	05-Aug-19 8:30:03 PM	141	96	73	62	A	
24	05-Aug-19 9:00:03 PM	152	107	85	77	A	
25	05-Aug-19 9:30:03 PM	133	81	55	61	E 2	Artifact/erratic oscillometric signal
26	05-Aug-19 9:35:52 PM	149	105	83	75	A	
27	05-Aug-19 10:30:03 PM	128	87	67	64	A	
28	05-Aug-19 11:30:03 PM	130	89	69	55	A	
29	06-Aug-19 12:30:03 AM	127	96	81	69	A	
30	06-Aug-19 1:30:03 AM	113	76	58	56	A	
31	06-Aug-19 2:30:03 AM	114	81	64	53	A	
32	06-Aug-19 3:30:03 AM	135	94	74	54	A	
33	06-Aug-19 4:30:03 AM	122	92	77	55	A	
34	06-Aug-19 5:18:42 AM	-	-	-	-	CH 112	
35	06-Aug-19 5:30:03 AM	128	95	78	48	A	
36	06-Aug-19 6:00:03 AM	-	-	-	-	E 87	Reading aborted—inflate time-out or airleak
37	06-Aug-19 6:04:21 AM	-	-	-	-	E 87	Reading aborted—inflate time-out or airleak
38	06-Aug-19 6:30:03 AM	126	92	75	56	A	
39	06-Aug-19 7:00:03 AM	130	92	73	61	A	
40	06-Aug-19 7:30:03 AM	131	100	84	61	A	
41	06-Aug-19 8:00:03 AM	131	95	77	60	A	
42	06-Aug-19 8:30:03 AM	128	96	80	55	A	
43	06-Aug-19 9:00:03 AM	133	98	81	53	A	
44	06-Aug-19 9:30:03 AM	149	106	84	55	A	
45	06-Aug-19 10:00:03 AM	152	105	82	54	A	
46	06-Aug-19 10:30:03 AM	151	108	87	46	A	
47	06-Aug-19 11:00:03 AM	142	107	90	52	A	
	06-Aug-19 11:15:56 AM	-	-	-	-	U 107	
	07-Aug-19 3:56:10 PM	-	-	-	-	U 110	

header

Patient ID	1234567	Weight	80 kg	Height	175 cm
MRN	6789	Address	Second street, 200, Big City		
Last Name	Doe	Order	555667		
First Name	John	Print Date	08-Mar-23		
Birth Date	01-Jan-70	Age	53	Test Date	05-Aug-19
Gender	Male	Recorder model	ABPM Device	S/N	ENG356

Physician Comments

Physician can wrote here some comments

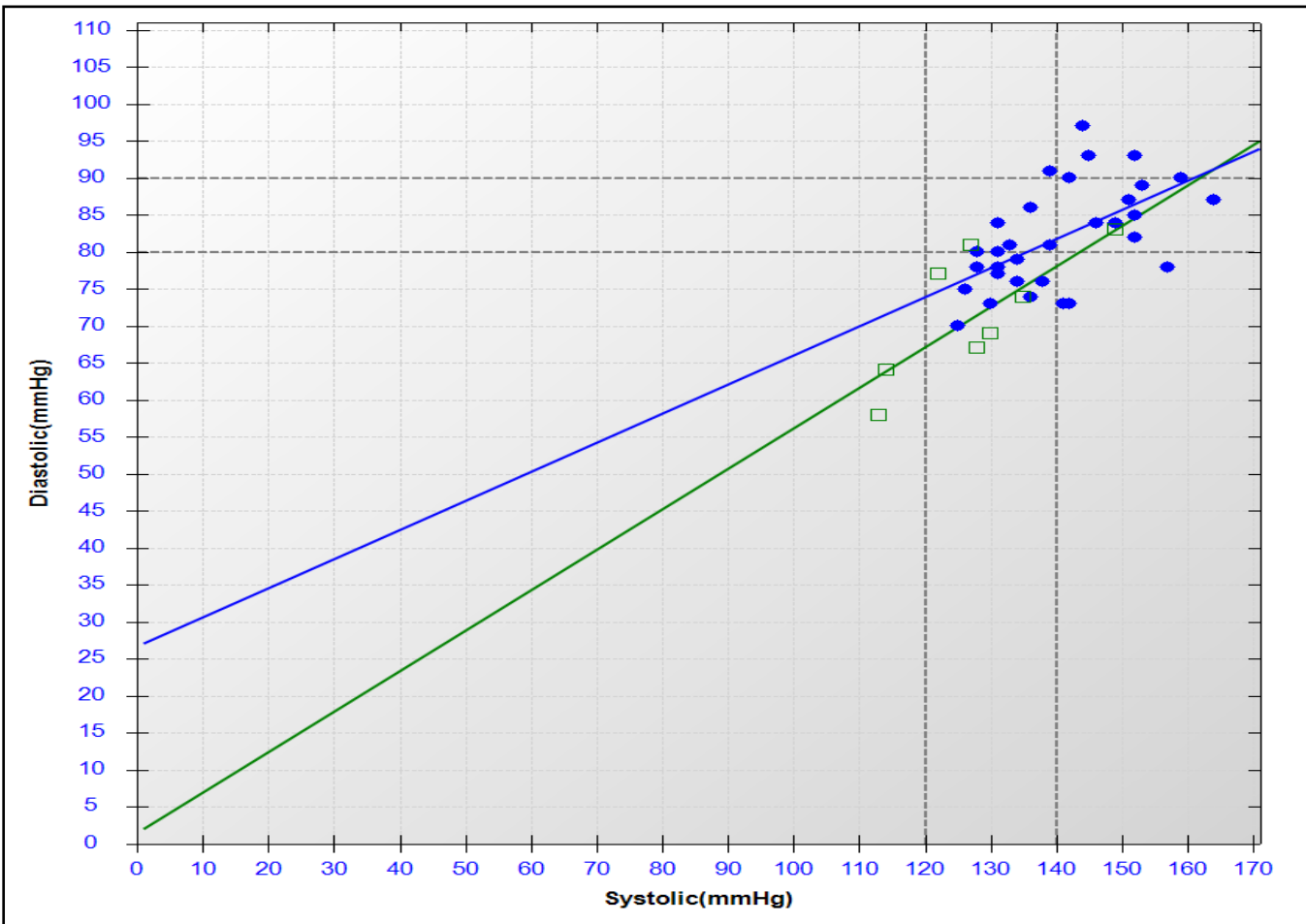

Awake Time ● Asleep Time □
 $DIA = 0.39SYS + 27.28$ — $DIA = 0.55SYS + 1.82$ —
 AASI 0.5 Correlation Coefficient $r = 0.71$ Average (SYS,DIA):(138,80)







header

Patient ID	1234567	Weight	80 kg	Height	175 cm
MRN	6789	Address	Second street, 200, Big City		
Last Name	Doe	Order	555667		
First Name	John	Print Date	08-Mar-23		
Birth Date	01-Jan-70	Age	53	Test Date	05-Aug-19
Gender	Male	Recorder model	ABPM Device	S/N	ENG356

Physician Comments

Physician can write here some comments

Values Above Limit

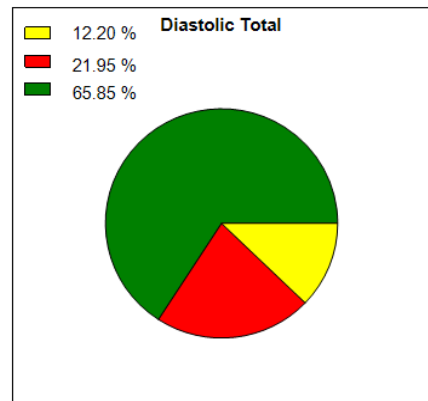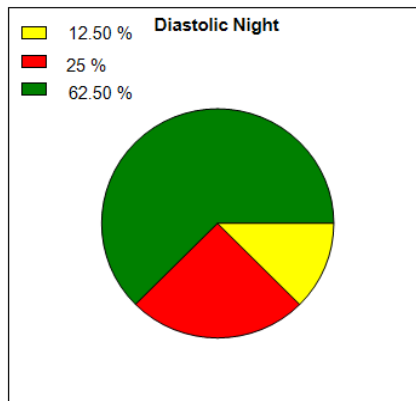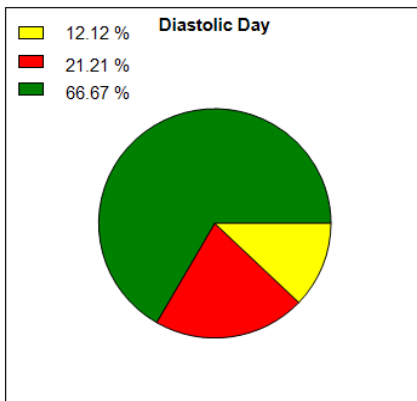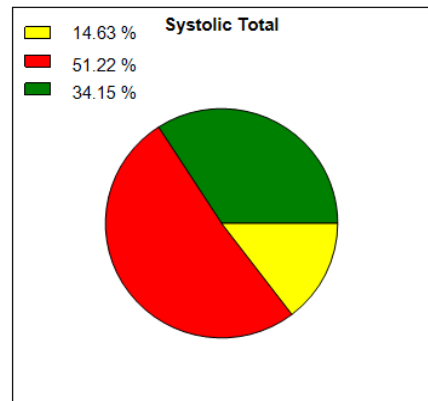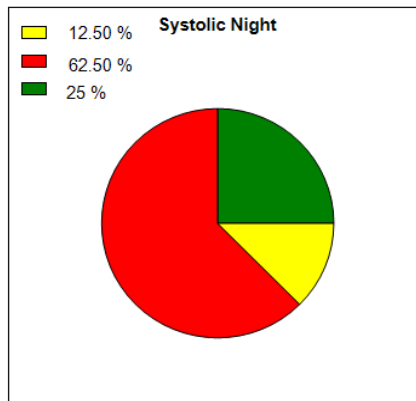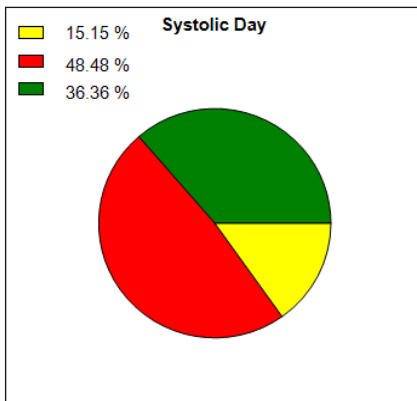

Acceptable Too High Normal

Patient ID	1234567	Weight	80	kg	Height	175	cm
MRN	6789	Address	Second street, 200, Big City				
Last Name	Doe	Order	555667				
First Name	John	Print Date	08-Mar-23				
Birth Date	01-Jan-70	Age	53	Test Date	05-Aug-19		
Gender	Male	Recorder model	ABPM Device	S/N	ENG356		

Summary

Period	Time	Samples	Valid Samples	Valid Samples (%)	BP Load Sys (%)	BP Load Dia (%)	Maximum SYS	Maximum DIA	Minimum SYS	Minimum DIA
Overall	05-Aug-19 10:37 06-Aug-19 12:00 (25:23)	48	41	85	54	22	164	97	113	58
Awake Period	05:18-21:30 (16:12)	40	33	82	48	21	164	97	125	70
Asleep Period	21:30-05:18 (07:48)	9	8	89	75	25	149	83	113	58
Dipping	SYS = 9.71% DIA = 13.01% <0% Inverted; <10% Non-Dipper; <20% Normal; >=20% Extreme									

Average Values

Measurements	Awake Period		Asleep Period		Overall	
	Value	Goal	Value	Goal	Value	Goal
SYS (mmHg)	141	<135	127	<120	138	<130
DIA (mmHg)	82	<85	72	<75	80	<80
MAP (mmHg)	102		90		99	
HR (BPM)	63		60		62	
PP (mmHg)	58	<60	55		58	