

# Burette Set

## 150mL In-Line Port

### DESCRIPTION

The m|devices Burette Set In-Line Port has a 150mL capacity vented chamber with a floating safety shut off valve. There is also a vented spike, roller regulator and a Swabable Needleless Access Site at one end and In-line Port with click clamp at the other.

### FEATURES & BENEFITS

- Vented perforator
- 150mL Burette
- Rigid roller regulator for accurate flow
- Swabable Needleless Access Site
- Floating safety shut off valve
- Click clamp
- In-line Port
- Peel open pouch
- Single use only


### PRODUCT DATA

- Burette capacity: 150mL
- Tubing length: 25cm
- Drip chamber: 20 drops /mL
- Priming volume: 1.8mL
- Sterile EO
- ARTG Number: 304860
- The product is MRI and CT suitable
- DEHP NO = Not detected (Result less than the Method Detection Limit EN 14372)
- Not made with Natural Rubber Latex
- Material: DEHP Free PVC tubing, PETG, PP & Silicone, POM
- Shelf Life: 3 Years

### ISO STANDARDS

- ISO 8536-5-2004 Infusion Equipment for Medical Use- Part 5: Burette infusion sets for single use, gravity feed
- ISO 8536-4-2010 Infusion Equipment for medical use - Part 4: infusion sets for single use, gravity feed

### PRODUCT DIAGRAM

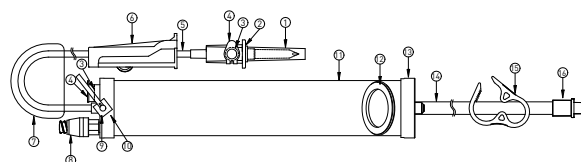

ITEM	DESCRIPTION
1	Spike Cap
2	Spike
3	Hydrophobic Membrane
4	Air-inlet
5	Tubing
6	Roller Clamp
7	Anti-kink Clip
8	Needle Free Valve
9	Hanger
10	Upper Cover
11	150mL Burette Chamber
12	Blue Float Valve
13	Lower Cover
14	Tubing
15	Robert Clamp
16	Adapter

Burette Set			STERILE
RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV040001	Burette Set	25cm	40