

Vented Double Spike Adaptor 30cm

DESCRIPTION

The Vented Double Spike Adaptor is designed with two vented perforators with roller regulators and a sealed septum connector. Made from medical grade non-toxic PVC (latex free).

FEATURES & BENEFITS

- 2 x Vented perforators
- 2 x Roller regulators
- Sealed septum connector
- Allows adaption from single to bag to double bag
- Peel open pouch


PRODUCT DATA

- Tubing diameter: 3mm
- Tubing length: 30cm
- Priming volume: 5.4mL
- Sterile EO
- The product is MRI and CT suitable
- ARTG Number: 304860
- Not Made with Natural Rubber Latex
- DEHP NO = Not detected (Result less than the Method Detection Limit EN 14372)
- Material: PVC
- Shelf Life: 3 Years

ISO STANDARDS

- ISO 8536-4 Infusion equipment for medical use - Part 4: Infusion sets for single use, gravity feed
- ISO 80369-7 Connector with 6% luer taper for intravascular or subcutaneous injection
- ISO 14971 Medical devices - Application of risk management to medical device

PRODUCT DIAGRAM

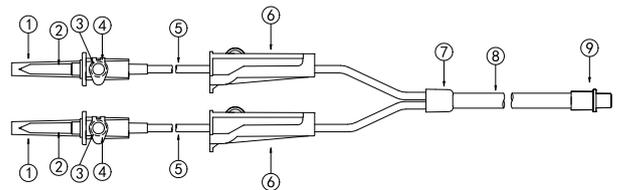

ITEM	DESCRIPTION
1	Cap
2	Perforator
3	Hydrophobic membrane
4	Air inlet cap
5	Tubing
6	Roller Clamp
7	Connector
8	Tubing
9	Adaptor

Vented Double Spike Adaptor			STERILE
RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV071001	Infusion/Transfusion Vented Double Spike Adaptor	30cm	50