

Mepilex® Border Post-Op

The flexible all-in-one post-op dressing


- Helps reduce the risk of surgical site infection (SSI) by minimising number of dressing changes and peri-wound skin damage¹⁻⁴
- Supports patient mobilisation by being flexible and gently adhering¹⁻⁴
- Reduce dressings related cost by minimising number of dressing changes¹⁻²

Supports longer wear times and low frequency of dressing changes¹⁻⁴

Shower-proof seal and bacteria barrier (microbes larger than 25 nm)

Ultra-absorbent material absorbs more blood than any other current leading dressing on the market¹¹⁻¹²

Transparent border allows inspection without removal

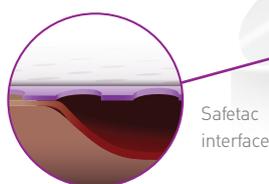

Clinically shown to minimise dressing-related skin damage and pain¹⁻⁴

Safetac interface adheres gently, maintaining skin integrity¹⁻⁴

Supports early patient mobilisation¹⁻⁴

Unique flex-cut pad stretches in all directions, conforming to the body even as it moves



Proven choice for a better outcome

Safetac*, pioneered by Mölnlycke Health Care, delivers above and beyond the ordinary. Proven to help optimise the wound healing journey and even prevent wounds, dressings with Safetac are the safe choice⁵ for patients and a champion for higher standards in wound care.

In fact, we have a wealth of evidence that supports the clinical and economic benefits of dressings with Safetac, including Mepilex®, Mepitel®, Mepiform® and Mepitac®. To date, these dressings have helped millions of patients worldwide⁵⁻⁷.

* A unique proprietary technology exclusive to Mölnlycke Health Care

How to use Mepilex® Border Post-Op



1. Open the sterile packaging and remove the dressing



2. Remove the middle part of the release film and apply the dressing on the right position



3. Remove the larger of the remaining films continuously while applying the dressing. Repeat for the smaller film and reposition if needed



4. Finalise the application by stroking the full dressing area for maximal adherence

How Mepilex Border Post-Op works

Mepilex Border Post-Op is an all-in-one post-op dressing that very effectively absorbs and retains surgical exudate. The Safetac interface minimise painful wound and peri-wound skin damage at dressing removal¹⁻⁴. The Safetac interface seals the wound edges, preventing the exudate to leak onto surrounding skin, minimising the risk of maceration⁸⁻¹². The flex-cut pad gives high flexibility and very good conformability over joints, such as knees or hips, promoting patient mobilisation¹⁻⁴.

Frequency of change

Mepilex Border Post-Op may be left in place for up to seven days depending on the condition on the wound and the surrounding skin, or as indicated by accepted clinical practice.

Benefits of Mepilex Border Post-Op

- Minimises skin damage, including blistering¹⁻⁴
- Very high absorption capacity leading to fewer dressing changes^{1-4, 10-12}
- High flexibility giving excellent comfort and conformability¹⁻³
- Can be lifted and adjusted without losing its adherent properties¹³
- Wide transparent borders for easy wound area inspection
- Bacteria and viral barrier (microbes >25nm)
- Leaves no residues low potential for skin irritation and allergy¹⁻⁴
- Shower proof¹⁻³

Mepilex Border Post-Op ordering information

Art. no	Size (cm)	Pcs/shelf	Pcs/transp.
496100	6x8	10	80
496200	9x10	10	70
496300	10x15	10	100
496400	10x20	10	120

Areas of use

Mepilex Border Post-Op is designed for exuding wounds. It is intended for acute wounds, such as surgical wounds, cuts and abrasions. It is optimised for post-op use and blood absorption. The design gives very high flexibility and makes it ideal to be used over joints such as hips and knees.

Precautions

In case of signs of clinical infection, consult a health care professional for adequate infection treatment.

Do not use on patients with known sensitivity to the dressing or its components.



Art. no	Size (cm)	Pcs/shelf	Pcs/transp.
496450	10x25	10	60
496600	10x30	10	40
496650	10x35	5	55

REFERENCES: 1. Zarghooni K. et al. Effect of a modern dressing compared to standard dressings on outcome after primary hip and knee arthroplasty: a prospective, non-randomised controlled study. E-poster, EWMA, London, 2015. 2. Bredow J. et al. Randomized clinical trial to evaluate performance of flexible self-adherent absorbent dressing coated with silicone layer after hip, knee or spinal surgery in comparison to standard wound dressing. Poster presentation at 5th Congress of the WUWHs, Florence, Italy, 2016. 3. Van Overschelde P. et al. A randomised controlled trial comparing two wound dressings used after elective hip and knee arthroplasty. Poster presentation at 5th Congress of the WUWHs, Florence, Italy, 2016. 4. Johansson C. et al. Preventing post-operative blisters following hip and knee arthroplasty. Wounds International, 2012. 5. Davies P. et al. Evidence review: the clinical benefits of Safetac technology in wound care. Journal of Wound Care, 2008. 6. Santamaria N. et al. Clinical effectiveness of a silicone foam dressing for the prevention of heel pressure ulcers in critically ill patients: Border II Trial. Journal of Wound Care, 2015. 7. Santamaria N. et al. An estimate of the potential budget impact of using prophylactic dressings to prevent hospital-acquired PUs in Australia. Journal of Wound Care, 2014. 8. Meaume S. et al. A study to compare a new self-adherent soft silicone dressing with a self-adherent polymer dressing in stage II pressure ulcers. Ostomy Wound Management, 2003. 9. Wiberg A.B. et al. Preventing maceration with a soft silicone dressing: in-vitro evaluations. Poster presented at the 3rd Congress of the WUWHs, Canada, 2008. 10. Feili F. et al. Blood absorption capacity of post-operative wound dressings. Poster presented at the 5th Congress of the WUWHs, Italy, 2016. 11. Feili F. et al. A laboratory evaluation of the fluid retention properties of post-operative absorbent dressings. Poster presented at the 5th Congress of the WUWHs, Italy, 2016. 12. Feili, F. et al. Fluid handling properties of post-operative wound dressings. Poster presentation at 5th Congress of the WUWHs, Florence, Italy, 2016. 13. Johnstone, A. et al. Innovations in the reduction of pressure ulceration and pain in critical care. Wounds UK, 2013.

Find out more at www.molnlycke.com.au

Mölnlycke Health Care, 12 Narabang Way, Level 4, Belrose, NSW, 2085. Phone 1800 005 231. New Zealand Orders & Enquiries 0800 005 231, www.molnlycke.co.nz. The Mölnlycke®, Mepilex®, Mepiform®, Mepitac®, Mepitel® and Safetac® trademarks, names and respective logos are registered globally to one or more of the Mölnlycke Health Care Group of Companies. © 2017 Mölnlycke Health Care. All rights reserved. AUWC0090

