

Dear valued customers,

The next generation of discs for your **cobas® b 101** is here.

Meet the new **cobas® click** test discs. They boast advanced performance and deliver reliable results. The discs are just a click away now.




cobas® click test discs for HbA1c, Lipids, and CRP

The most important innovations at a glance.

All Discs:

- ✓ **Improved disc's sample inlet.**
You can now place the discs on the blood drop from all sides:

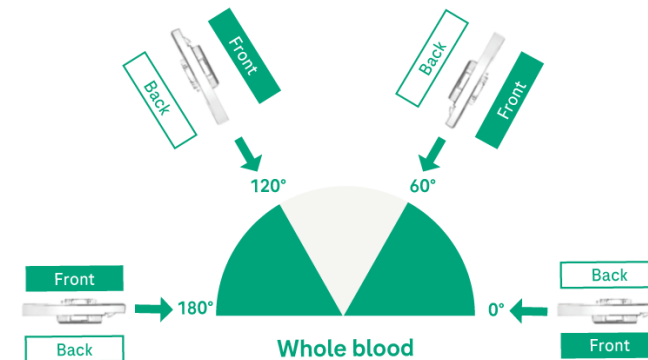

- ✓ **Just one click to close**
- ✓ **Optimised package size:**
Now with **20 test discs** instead of 10 per package

Lipid Discs:

- ✓ **New failsafe method for triglycerides.**
Detects interfering substances such as hand cream residues at the blood drawing facility and adjusts the TG value. A measurement error is displayed if the error is too severe and cannot be rectified. In addition, the reduced sample size is now 17µl instead of 19µL¹.

HbA1c Discs:

- ✓ **Optimised HbA1c performance.**
Improved denaturation in HbA1c test to minimise Hb variant interferences².

CRP Discs:

- ✓ **No adjustment to performance or sample time.**

¹ - Roche Diagnostics International Ltd. cobas click Lipid Panel Test Method Sheet.

² - Roche Diagnostics International Ltd. cobas HbA1c Test Gen. 2 Method Sheet.



You will need to update the software to version 3.0.0 to use the new **cobas® click** discs on your **cobas® b 101**. If you have not yet installed this update, please do so before the upgrade. Please note that the previous **cobas® b 101** discs will be discontinued at the beginning of 2026. We kindly request you to switch to the new **cobas® click** discs in due time.

COMMENCING JUNE 2026

From June 2026 and onwards, please order the new discs 2.0 with the following item numbers:

09200509190 **cobas®** HbA1c Test Gen. 2 instead of 08038694190
09200495190 **cobas®** Lipid Panel Gen. 2 instead of 06380115190
09200487190 **cobas®** CRP Test Gen. 2 instead of 08024669190

Moreover, they come with the new controls.

These are the same and readily available:

09355537190 **cobas®** HbA1c Control Gen. 2 instead of 06380204190
09355545190 **cobas®** Lipid Control Gen. 2 instead of 06380182190
09355553190 **cobas®** CRP Control Gen. 2 instead of 08024723190



cobas® click System

In addition to the new **cobas® click** test discs, we are also pleased to introduce the new **cobas® click** instrument, the successor to the **cobas® b 101**.

- ✓ Features an integrated barcode scanner
- ✓ Enhanced connectivity including WLAN
- ✓ Capacitive touch screen for more ease of use
- ✓ Larger memory capacity

Should you have any questions about **cobas® click** or the software update for **cobas® b 101**, please contact our customer support by email at australia.dia_cardiology@roche.com or phone **1800 645 619**.



Learn more about the new **cobas® click** and the **cobas® click** test discs. Both are now available.

COBAS and COBAS CLICK are trademarks of Roche.
© 2026 Roche

Roche Diagnostics Australia Pty Limited
ABN 29 003 001 205
2 Julius Avenue, North Ryde NSW 2113
Phone: 02 9860 2222
MC-AU-04418

THIS PRODUCT IS NOT AVAILABLE FOR
PURCHASE BY THE GENERAL PUBLIC.
ALWAYS READ THE LABEL AND FOLLOW
THE DIRECTIONS FOR USE.