

Welch Allyn Green Series™ Diagnostic Desk System

Now, everything you need for comprehensive patient exams is within reach—for a simplified workflow and faster patient throughput.

Part number 77535-C1 shown—includes the Green Series Diagnostic Desk System, desk aneroid w/adult cuff, 3.5 V Ni-Cad handle with diagnostic otoscope and Braun ThermoScan® PRO 4000 thermometer.



WelchAllyn®



OTOSCOPE



OPHTHALMOSCOPE




THERMOMETER

Welch Allyn Green Series Diagnostic Desk System provides you with
all of the tools you need in one compact set,
 helping to reduce exam time. And you can choose the specific tools that work
 best for your facility, maximizing your investment. Plus, the Green Series™ name
 promises you improved energy savings, to improve your bottom line.

3.5 V Diagnostic Otoscope*

- Halogen HPX® lamp provides light output for true tissue color and long-lasting performance
- Fiber optics project cool light to the distal tip with no reflections or obstructions for superior views of the tympanic membrane
- Wide-angle viewing lens allows for instrumentation under magnification
- Sealed system for pneumatic otoscopy and accurate diagnosis of otitis media

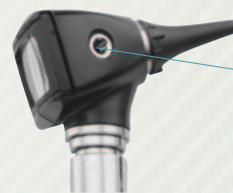

Insufflation Port

MacroView™ Otoscope*

- A nearly complete view of the tympanic membrane with approximately twice the field of view and 30% greater magnification than a traditional otoscope
- Better clarity and definition of landmarks
- Provides the ability to adjust focus for variable ear canal length and farsighted eyes



Tip Grip: Ensures ear specula are fastened securely and easily

Insufflation Port for Pneumatic Otoscopy

Default Focus: Optimal setting for most ear examinations

Adjustable Focus: Ability to zoom in or out to fine-tune the view

Welch Allyn 3.5 V Coaxial Ophthalmoscope*

Welch Allyn 3.5 V coaxial ophthalmoscope allows for easy entry into the eye, a larger field of view and reduced glare compared to standard ophthalmoscopes.

- Halogen HPX lamp provides light for true tissue color and long-lasting performance
- Produces a shadow-free spot for easier entry into undilated pupils



Triple-sealed Ophthalmoscope

Coaxial Optical System

Cross-linear Polarizing Filter

Aperture Selection Dial

Lens Selection Disc

PanOptic™ Ophthalmoscope*

- Fast, easy entry into undilated pupils
- 5X larger view of the fundus vs. standard ophthalmoscopes in an undilated eye
- 25° field of view vs. the standard 5° field of view



Axial PointSource™ Optics

Patient Eye Cup

Patented Glare Extinguishment System

Dynamic Focusing Wheel

Soft Grip Handle

Easy-to-use Aperture Dial

Each desk system comes with a Welch Allyn Braun ThermoScan PRO 4000 Ear Thermometer and a Desk Aneroid with an Adult FlexiPort® Blood Pressure Cuff



Braun ThermoScan® PRO 4000

- Provides the speed and accuracy you need for all types of patients through a preheated probe tip and sensor
- ExacTemp™ technology provides feedback to the clinician to ensure proper positioning and accuracy
- Probe cover detection system ensures precise infrared heat transfer for accurate temperatures



Desk Aneroid

- Laser-engraved dial face provides pinpoint accuracy
- Lifetime calibration warranty

Specifications: Welch Allyn Green Series™ Diagnostic Desk System

Environment	Operating: Temperature Range: 59°F (15°C)–104°F (40°C) Humidity Range: 15%–95% Atmospheric Pressure Range: 700hPa–1060hPa
	Transport/Storage: Temperature Range: 4°F (-20°C)–120°F (50°C) Humidity Range: 15%–95% Atmospheric Pressure Range: 700hPa–1060hPa
Electrical	Power Supply: AC Input Power: 100–240 V, 50/60 Hz, 0.18 A Max DC Output Power: 5 V, 1 A
	Charging Device: DC Input Power: 5 V, 1 A DC Output Power: 5 V, 0.5 A
Warranty	Green Series Diagnostic Desk System (One Year) Desk Aneroid with Adult Blood Pressure Cuff (Lifetime Calibration) 3.5 V Ni-Cad Handle (Two Years) (2) 3.5 V Lithium Ion Handles (One Year) Diagnostic Otoscope Head (One Year) 3.5 V Coaxial Ophthalmoscope Head (One Year) MacroView™ Otoscope Head w/ Throat Illuminator (One Year) PanOptic Ophthalmoscope Head (One Year) Braun ThermoScan® PRO 4000 Infrared Ear Thermometer (Three Years Limited Warranty)

Ordering Information

77535-C1 77535-14 77535-17	Green Series™ Diagnostic Desk System Desk Aneroid w/ Adult Blood Pressure Cuff 3.5 V Ni-Cad Handle Diagnostic Otoscope Head Braun ThermoScan® PRO 4000 Infrared Ear Thermometer
77535-C2 77535-26 77535-27	Green Series Diagnostic Desk System Desk Aneroid w/ Adult Blood Pressure Cuff 3.5 V Ni-Cad Handle MacroView™ Otoscope Head w/ Throat Illuminator Braun ThermoScan® PRO 4000 Infrared Ear Thermometer
77535-C3 77535-34 77535-37	Green Series Diagnostic Desk System Desk Aneroid w/ Adult Blood Pressure Cuff 3.5 V Lithium Ion Handle Diagnostic Otoscope Head Braun ThermoScan® PRO 4000 Infrared Ear Thermometer
77535-C4 77535-46 77535-47	Green Series Diagnostic Desk System Desk Aneroid w/ Adult Blood Pressure Cuff 3.5 V Lithium Ion Handle MacroView™ Otoscope Head w/ Throat Illuminator Braun ThermoScan® PRO 4000 Infrared Ear Thermometer
77535-C5 77535-54 77535-56	Green Series Diagnostic Desk System Desk Aneroid w/ Adult Blood Pressure Cuff (2) 3.5 V Ni-Cad Handles Diagnostic Otoscope Head Coaxial Ophthalmoscope Head Braun ThermoScan® PRO 4000 Infrared Ear Thermometer
77535-66 (Australia only)	Green Series Diagnostic Desk System Desk Aneroid w/ Adult Blood Pressure Cuff (2) 3.5 V Lithium Ion Handles MacroView™ Otoscope Head w/ Throat Illuminator 3.5 V Coaxial Ophthalmoscope Head Braun ThermoScan® PRO 4000 Infrared Ear Thermometer
77535-76 (Australia only)	Green Series Diagnostic Desk System w/ Chargers Only
77535-86 (Australia only)	Green Series Diagnostic Desk System Desk Aneroid with Adult Blood Pressure Cuff (2) 3.5 V Lithium Ion Handles MacroView™ Otoscope Head w/ Throat Illuminator PanOptic Ophthalmoscope Head Braun ThermoScan® PRO 4000 Infrared Ear Thermometer

Australia: +61 2 9638 3000
China: +86 512 68086078
Japan: +813 3219 0071

Malaysia: +603 788 43341
Singapore: +65 6419 8100
South Africa: +27 11 777 7555

WelchAllyn®